




TrimTM

The first truly accessible
professional aligner trimming platform.

The Final Manual Step in Aligner Manufacturing — Now Fully Automated and Accessible to Everyone.

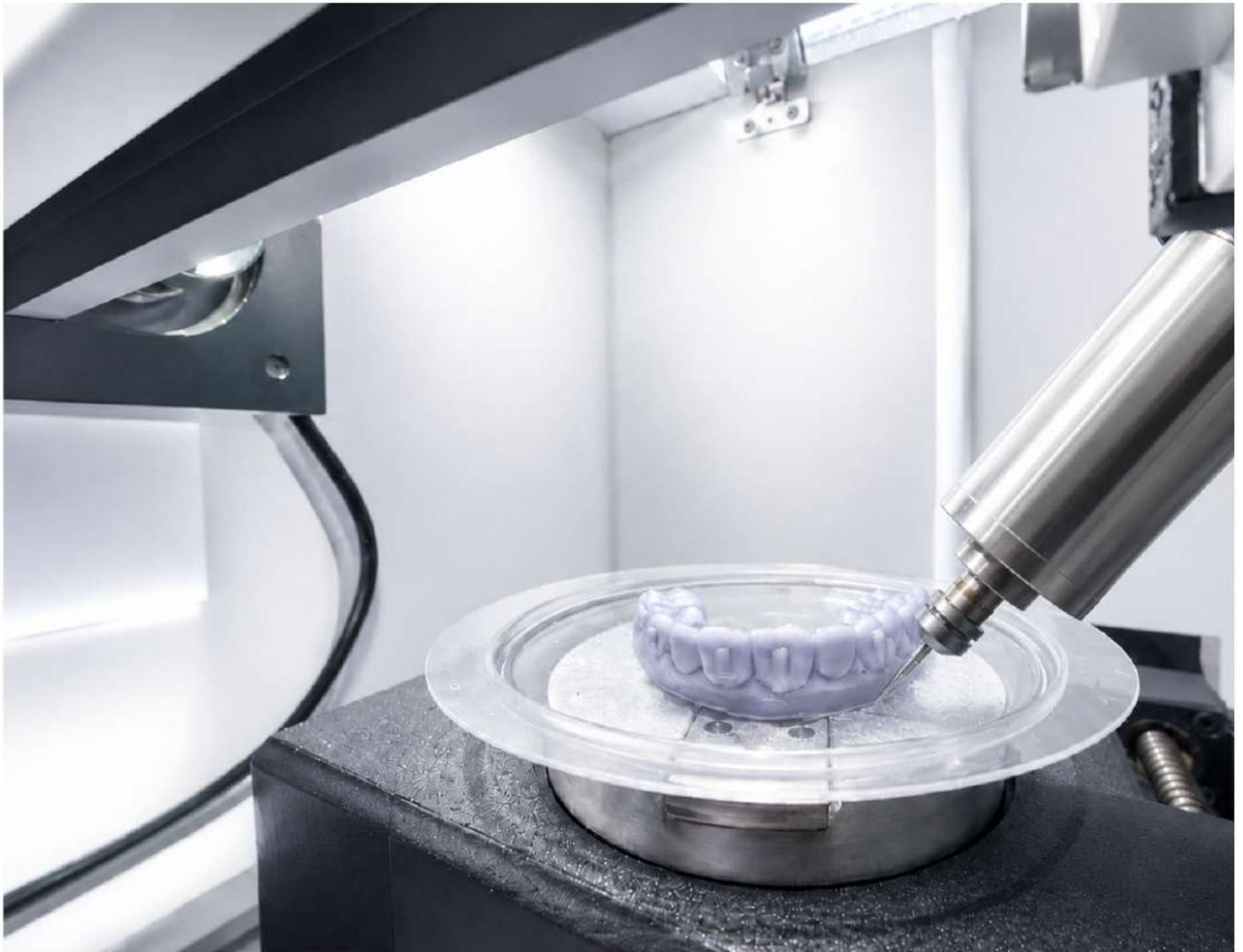

Traditional robotic trimming systems were often expensive, closed, and restricted to specific software ecosystems. Our approach is different. The Trim Automatic Aligner Trimmer is built as a completely open platform, allowing labs and clinics to integrate automated trimming into their existing workflow without changing their preferred aligner software. If your software can export DPA or aligner shell files, the system is ready to work !

How Trim[™] made this possible?

Proprietary software


Trim is powered by a fully in-house developed CAM engine through its dedicated software ecosystem — TrimBase and TrimCAM — eliminating dependence on third-party milling software while enabling seamless integration between hardware, trimline generation, and precision toolpath optimization.

Inhouse Technology

Because Trim is built on a completely in-house technology stack, we eliminate the high licensing and integration costs associated with third-party systems. From machine design to electronics, every part of the workflow is developed internally — allowing us to optimize performance, reduce dependency costs, and deliver advanced automation at a significantly more accessible price point.



Proprietary fixture software for robotic trimming

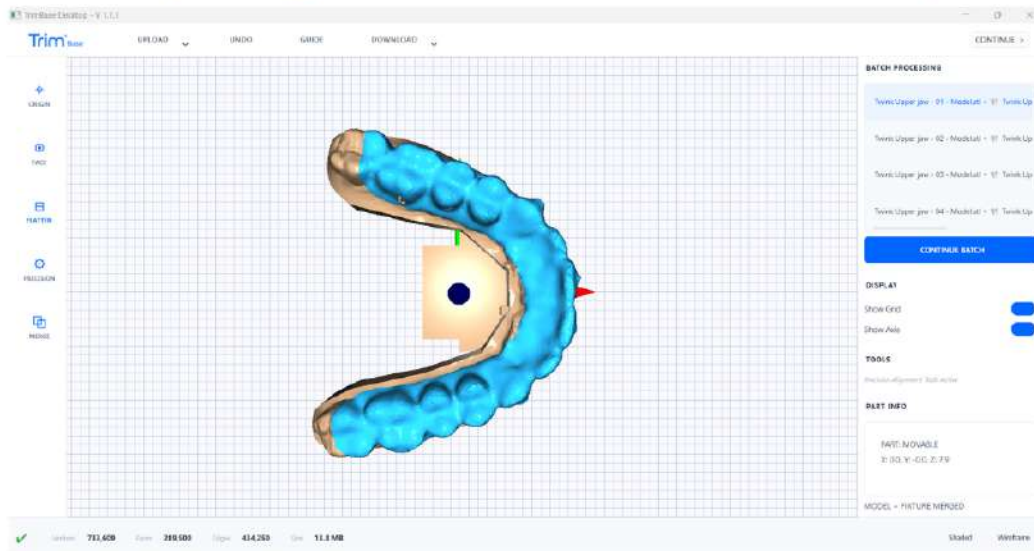


TrimTMBase

TrimBase is an advanced software solution designed for automated aligner model trimming by efficiently adding customizable fixtures to dental models in batch built for seamless integration with TRIM automated trimming machines.

Proprietary CAM engine for robotic trimming

TrimTMCAM

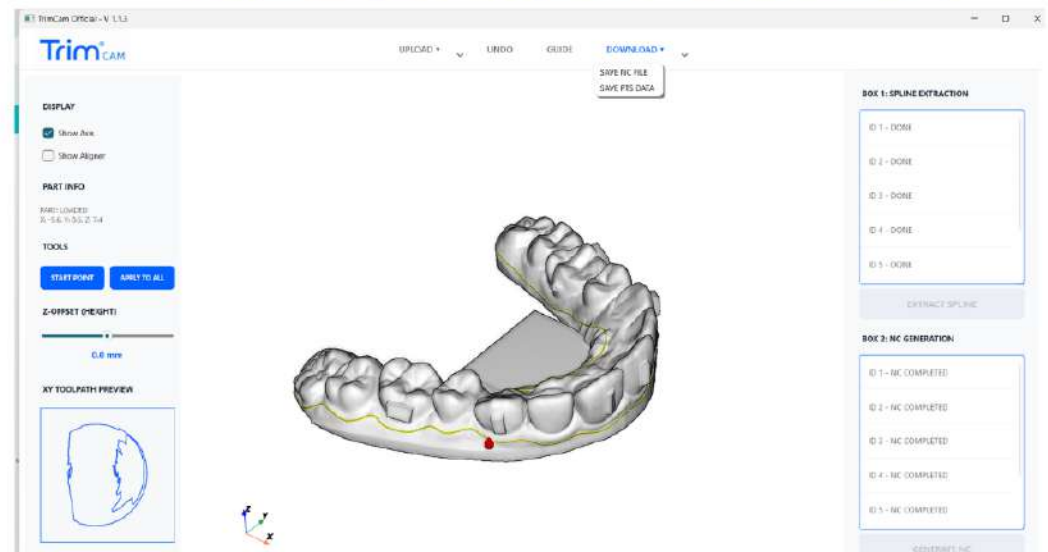


TrimCAM is an intelligent software that derive trimlines directly from dental models and aligners, allowing users to generate accurate trimming boundaries with minimal manual intervention. In addition, TrimBase features a powerful integrated CAM engine capable of generating precise toolpaths for 4-axis trimming machines, ensuring reliable automated machining performance.

Easy release fixture

Our system is equipped with an easy-release fixture system that allows models to be attached and removed with a simple press — eliminating the need for screws or manual fastening.

This enables faster loading and unloading, improves accessibility during operation, and significantly reduces overall trimming time.

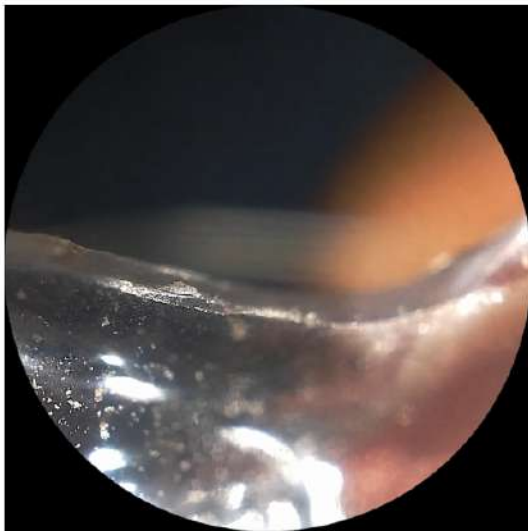

Brushless motor spindle

Driven by a high-speed 36,000 RPM Spindle, Trim achieves extremely refined edge finishing with remarkable consistency — reducing post-processing and requiring only minimal polishing before delivery.

Compact and Tabletop

Engineered with a compact tabletop form factor of just 410mm × 365mm × 470mm (L × B × H), Trim integrates effortlessly into modern dental labs and clinics. Direct PC connectivity enables seamless file transfer, live monitoring, and an efficient digital workflow experience.



Specifications

Specification	Details
Machine Dimensions	410mm × 365mm × 470mm (L × B × H)
Machine Weight	20 kg
Spindle	36,000 RPM Spindle
Release Mechanism	Quick-Release 3-Point Fixture System
Average Trimming Time	1 Minute 15 Seconds
Connectivity	USB Type-A via PC
Power Supply	110V / 220V Compatible
Dust Collection	External Vacuum Outlet Provision
Tool Holding System	Manual Twist-Lock Chuck



TrimTM

